

Human Osteopontin ELISA Instructions

CAT: AEH0057

CONTENT

No.	Content	CAT. No	Volume
1	CP (Coated Plate)	EH0057CP	96 wells
2	S (Standard)	EH0057S1	2 vials
3	SD (Sample Diluent)	ESD06	15 ml/bottle
4	DA (Detect Antibody)	EH0057DA	6 ml/bottle
5	SH (Streptavidin-HRP)	ESH01	12 ml/bottle
6	AB (Assay Buffer 1×)	EAB01	12 ml/bottle
7	TS (TMB Substrate)	ETS01	12 ml/bottle
8	SS (Stop Solution)	ESS01	12 ml/bottle
9	WB (Wash Buffer 10×)	EWB01	50 ml/bottle
10	SF (Sealer Film)	ESF01	6 pieces

NOTE: After the kit is opened, the stabilization period of each content is 30 days, so please use the kit within 30 days after opening.


REAGENT PREPARATION

Washing Buffer (1×) Preparation: Pour entire contents (50 ml) of the **Washing Buffer Concentrate** (10×) into a clean 500 ml graduated cylinder. Bring to final volume of 500 ml with glass-distilled or deionized water. Transfer to a clean wash bottle and store at 2 to 25°C.

Standard Curve Preparation: Reconstitute Human OPN Standard by addition of distilled water as S1. Reconstitution volume is stated on the label of the standard vial. Swirl or mix gently to insure complete and homogeneous solubilization (concentration of reconstituted standard = 20 ng/ml).

Allow the standard to reconstitute for 10-30 minutes. Mix well prior to making dilutions.

Pipette 200 µl of Sample Diluent into each tube. Use the high standard to produce a 1:1 dilution series. Mix each tube thoroughly before the next transfer. Sample Diluent serves as the zero standard (0 ng/ml).

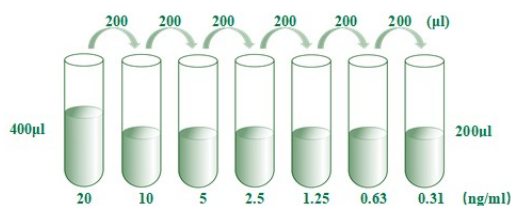


ASSAY PROCEDURE

Bring all reagents and samples to room temperature before use.

1. Prepare all reagents and working standards as directed in the previous sections.
2. Remove excess **CP** (Coated Plate) strips from the plate frame, return them to the foil pouch and reseal.
3. Add 50 µl of **AB** (Assay Buffer) to each well.
4. Add 50 µl or 10 µl of **Standard or sample** per well. Ensure reagent addition is uninterrupted and completed within 15 minutes.
5. Add 50 µl of **DA** (Detect Antibody) to each well.
6. Cover with an **SF** (Sealer Film). Incubate at room temperature (18 to 25°C) for 1 hours on a microplate **shaker** set at 500 rpm.
7. Aspirate each well and **wash**, repeating the process four times. Wash by filling each well with WB (Washing Buffer 300 µl). Complete removal of liquid at each step is essential to good performance. After the last wash, remove any remaining WB (Washing Buffer) by aspirating or decanting. Invert the plate and **blot** it against clean paper towels.
8. Add 100 µl of **SH** (Streptavidin-HRP) to each well.
9. Cover with a new **SF** (Sealer Film). Incubate at room temperature (18 to 25°C) for 30 min on a microplate **shaker** set at 500 rpm.
10. Repeat aspiration/**wash** as in step 7.
11. Add 100 µl of **TS** (TMB Substrate) to each well. Incubate for 5-30 minutes at room temperature.
12. Add 100 µl of **SS** (Stop Solution) to each well.
13. Determine the optical density within 30 minutes, using microplate **reader** set to 450 nm corrected with 570 nm or 630 nm.

TYPICAL DATA

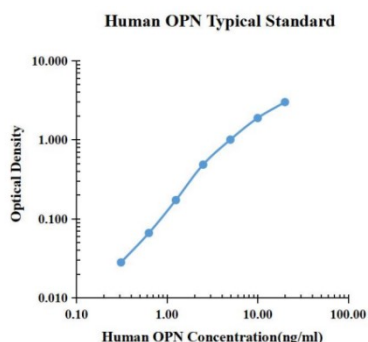

ng/ml	O.D.		Average	Corrected
0.00	0.0246	0.0244	0.0245	
0.31	0.0543	0.0502	0.0523	0.0278
0.63	0.0937	0.0866	0.0902	0.0657
1.25	0.2017	0.1903	0.1960	0.1715
2.50	0.5132	0.5027	0.5080	0.4835
5.00	1.0570	0.9935	1.0253	1.0008
10.00	1.9270	1.8780	1.9025	1.8780
20.00	3.0580	2.9620	3.0100	2.9855

SENSITIVITY

The minimum detectable dose (MDD) of human OPN is typically less than 0.04 ng/ml (50 µl of sample volume) or 0.12 ng/ml (10 µl of sample volume).

The MDD was determined by adding two standard deviations to the mean optical density value of ten zero standard replicates and calculating the corresponding concentration.

PRECISION

- Intra-assay Precision (Precision within an assay) Three samples of known concentration were tested twenty times on one plate to assess intra-assay precision.
- Inter-assay Precision (Precision between assays).

	Intra-assay Precision			Inter-assay Precision		
	S1	S2	S3	S1	S2	S3
Sample Number	22	22	22	6	6	6
Average (pg/ml)	0.42	2.09	6.22	0.51	2.10	6.03
Standard deviation	0.02	0.09	0.48	0.02	0.10	0.37
Coefficient of variation (%)	4.9	4.5	7.7	4.3	4.7	6.1

RECOVERY

The spike recovery was evaluated by spiking 3 levels of human OPN into health human serum sample. The un-spiked serum was used as blank in this experiment.

The recovery ranged from 83% to 117% with an overall mean recovery of 108%.

LINEARITY

To assess the linearity of the assay, five samples were spiked with high concentration of OPN in human serum and diluted with Sample Diluent to produce samples with values within the dynamic range of the assay.

The linearity ranged from 84% to 95% with an overall mean recovery of 89%.

SAMPLE VALUES

Serum/Plasma - Thirty samples from apparently healthy volunteers were evaluated for the presence of human OPN in this assay. No medical histories were available for the donors.

Sample Matrix	Sample Evaluated	Range (pg/ml)	Detectable %	Mean of Detectable (pg/ml)
Serum	30	1.10-6.74	100	3.97

n.d. = non-detectable. Samples measured below the sensitivity are considered to be non-detectable.

DISCLAIMER AND VERSION

For research use only. Not for use in clinical diagnostic procedures.

Copyright 2024 AREX Biosciences Ltd. Version.240820