

Human Osteoprotegerin/TNFRSF11B ELISA Instructions

CAT: AEH0070

CONTENT

No.	Content	CAT. No	Volume
1	CP (Coated Plate)	EH0070CP	96 wells
2	S (Standard)	EH0070S,S1~S7,S0	9 vials
3	DA (Detect Antibody)	EH0070DA	6 ml/bottle
4	SH (Streptavidin-HRP)	ESH01	12 ml/bottle
5	AB (Assay Buffer 1×)	EAB01	12 ml/bottle
6	TS (TMB Substrate)	ETS01	12 ml/bottle
7	SS (Stop Solution)	ESS01	12 ml/bottle
8	WB (Wash Buffer 10×)	EWB01	50 ml/bottle
9	SF (Sealer Film)	ESF01	6 pieces

NOTE: After the kit is opened, the stabilization period of each content is 30 days, so please use the kit within 30 days after opening.


REAGENT PREPARATION

Washing Buffer (1×) Preparation: Pour entire contents (50 ml) of the **Washing Buffer Concentrate** (10×) into a clean 500 ml graduated cylinder. Bring to final volume of 500 ml with glass-distilled or deionized water. Transfer to a clean wash bottle and store at 2 to 25°C.

Standard Curve Preparation: S1 to S7 and S0 is ready to use for serum and plasma.

Other sample type, prepare the standard curve with whatever buffer (SPB, Sample Prepared Buffer) is used to prepare the sample, such as cell culture supernatant, tissue grinding liquid, cell lysate, etc. Urine sample use AB (Assay Buffer) prepare standard curve.

The human Osteoprotegerin Standard EH0070S 50000 pg/ml 30 μ l + 270 μ l SPB serves as the high standard (5000 pg/ml). Pipette 200 μ l of SPB into each tube. Use the high standard to produce a 1:2 dilution series. Mix each tube thoroughly before the next transfer. SPB serves as the zero standard (0 pg/ml).

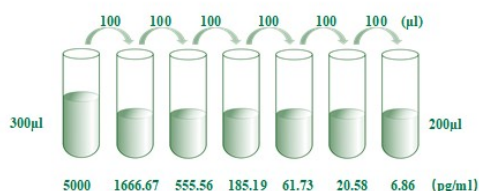


ASSAY PROCEDURE

Bring all reagents and samples to room temperature before use.

1. Prepare all reagents and working standards as directed in the previous sections.
2. Remove excess **CP** (Coated Plate) strips from the plate frame, return them to the foil pouch and reseal.
3. Add 50 μ l of **AB** (Assay Buffer) to each well.
4. Add 50 μ l or 10 μ l of **Standard or sample** per well. Ensure reagent addition is uninterrupted and completed within 15 minutes.
5. Add 50 μ l of **DA** (Detect Antibody) to each well.
6. Cover with an **SF** (Sealer Film). Incubate at room temperature (18 to 25°C) for 1 hours on a microplate **shaker** set at 500 rpm.
7. Aspirate each well and **wash**, repeating the process four times. Wash by filling each well with WB (Washing Buffer 300 μ l). Complete removal of liquid at each step is essential to good performance. After the last wash, remove any remaining WB (Washing Buffer) by aspirating or decanting. Invert the plate and **blot** it against clean paper towels.
8. Add 100 μ l of **SH** (Streptavidin-HRP) to each well.
9. Cover with a new **SF** (Sealer Film). Incubate at room temperature (18 to 25°C) for 30 min on a microplate **shaker** set at 500 rpm.
10. Repeat aspiration/**wash** as in step 7.
11. Add 100 μ l of **TS** (TMB Substrate) to each well. Incubate for 5-30 minutes at room temperature.
12. Add 100 μ l of **SS** (Stop Solution) to each well.
13. Determine the optical density within 30 minutes, using microplate **reader** set to 450 nm corrected with 570 nm or 630 nm.

TYPICAL DATA

Human Osteoprotegerin Typical Standard

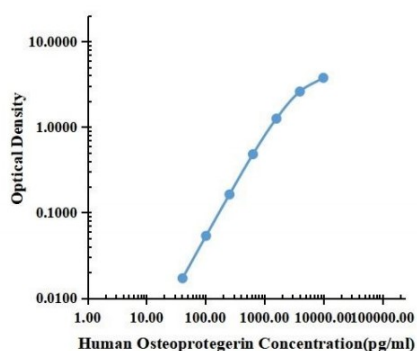

pg/ml	O.D.		Average	Corrected
0.00	0.0149	0.0139	0.0144	
6.86	0.0332	0.0309	0.0321	0.0171
20.58	0.0686	0.0681	0.0684	0.0534
61.73	0.1782	0.1770	0.1776	0.1626
185.19	0.4947	0.4959	0.4953	0.4803
555.56	1.2700	1.2650	1.2675	1.2525
1666.67	2.6800	2.5420	2.6110	2.5960
5000.00	3.8430	3.6970	3.7700	3.7550

SENSITIVITY

The minimum detectable dose (MDD) of human Osteoprotegerin is typically less than 0.37 pg/ml (50 μ l of sample volume) or 2.53 pg/ml (10 μ l of sample volume).

The MDD was determined by adding two standard deviations to the mean optical density value of ten zero standard replicates and calculating the corresponding concentration.

PRECISION

- Intra-assay Precision (Precision within an assay) Three samples of known concentration were tested twenty times on one plate to assess intra-assay precision.
- Inter-assay Precision (Precision between assays).

	Intra-assay Precision			Inter-assay Precision		
	S1	S2	S3	S1	S2	S3
Sample Number	22	22	22	6	6	6
Average (pg/ml)	194.3	749.0	1558.3	209.4	735.1	1459.6
Standard deviation	12.2	35.0	105.3	9.6	24.8	29.6
Coefficient of variation (%)	6.3	4.7	6.8	4.6	3.4	2.0

RECOVERY

The spike recovery was evaluated by spiking 3 levels of human Osteoprotegerin into health human serum sample. The un-spiked serum was used as blank in this experiment.

The recovery ranged from 65% to 98% with an overall mean recovery of 95%.

LINEARITY

To assess the linearity of the assay, five samples were spiked with high concentration of Osteoprotegerin in human serum and diluted with Sample Diluent to produce samples with values within the dynamic range of the assay.

The linearity ranged from 69% to 86% with an overall mean recovery of 72%.

SAMPLE VALUES

Serum/Plasma - Thirty samples from apparently healthy volunteers were evaluated for the presence of Osteoprotegerin in this assay. No medical histories were available for the donors.

Sample Matrix	Sample Evaluated	Range (pg/ml)	Detectable %	Mean of Detectable (pg/ml)
Serum	30	341.28-586.42	100	461.05

n.d. = non-detectable. Samples measured below the sensitivity are considered to be non-detectable.

DISCLAIMER AND VERSION

For research use only. Not for use in clinical diagnostic procedures.

Copyright 2024 AREX Biosciences Ltd. Version.240820