

## HIV-1 p24 Instructions

CAT: AEH0078

### CONTENT

No.	Content	CAT. No	Volume
1	CP (Coated Plate)	EH0078CP	96 wells
2	S (Standard)	EH0078S1,S1-S7,S0	9 vials
3	SD (Sample Diluent)	ESD01	12 ml/bottle
4	AB (Assay Buffer 1×)	EAB01	6 ml/bottle
5	DA (Detect Antibody )	EH0078DA	6 ml/bottle
6	SH (Streptavidin-HRP)	ESH01	12 ml/bottle
7	TS (TMB Substrate)	ETS01	12 ml/bottle
8	SS (Stop Solution)	ESS01	12 ml/bottle
9	WB (Wash Buffer 10×)	EWB01	50 ml/bottle
10	SF (Sealer Film)	ESF01	6 pieces

**NOTE:** After the kit is opened, the stabilization period of each content is 30 days, so please use the kit within 30 days after opening.


### REAGENT PREPARATION

**Washing Buffer (1×) Preparation:** Pour entire contents (50 ml) of the **Washing Buffer Concentrate** (10×) into a clean 500 ml graduated cylinder. Bring to final volume of 500 ml with glass-distilled or deionized water. Transfer to a clean wash bottle and store at 2 to 25°C.

**Standard Curve Preparation:** S1 to S7 and S0 is ready to use for serum, plasma and cell culture supernatant.

Other sample type, prepare the standard curve with whatever buffer (SPB, Sample Prepared Buffer) is used to prepare the sample, such as tissue grinding liquid, cell lysate, etc. Urine sample use AB (Assay Buffer) prepare standard curve.

The HIV-1 p24 Standard EH0078S 2000 pg/ml 30 µl + 270 µl SPB serves as the high standard (5000 pg/ml). Pipette 150 µl of SPB into each tube. Use the high standard to produce a 1:1 dilution series. Mix each tube thoroughly before the next transfer. SPB serves as the zero standard (0 pg/ml).

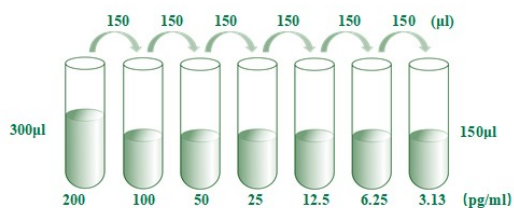


### ASSAY PROCEDURE

Bring all reagents and samples to room temperature before use.

1. Prepare all reagents and working standards as directed in the previous sections.
2. Remove excess **CP** (Coated Plate) strips from the plate frame, return them to the foil pouch and reseal.
3. Add 50 µl of **AB** (Assay Buffer) to each well.
4. Add 50 µl or 10 µl of **Standard or sample** per well. Ensure reagent addition is uninterrupted and completed within 15 minutes.
5. Add 50 µl of **DA** (Detect Antibody) to each well.
6. Cover with an **SF** (Sealer Film). Incubate at room temperature (18 to 25°C) for 1 hours on a microplate **shaker** set at 500 rpm.
7. Aspirate each well and **wash**, repeating the process four times. Wash by filling each well with WB (Washing Buffer 300 µl ). Complete removal of liquid at each step is essential to good performance. After the last wash, remove any remaining WB (Washing Buffer) by aspirating or decanting. Invert the plate and **blot** it against clean paper towels.
8. Add 100 µl of **SH** (Streptavidin-HRP) to each well.
9. Cover with a new **SF** (Sealer Film). Incubate at room temperature (18 to 25°C) for 30 min on a microplate **shaker** set at 500 rpm.
10. Repeat aspiration/**wash** as in step 7.
11. Add 100 µl of **TS** (TMB Substrate) to each well. Incubate for 5-30 minutes at room temperature.
12. Add 100 µl of **SS** (Stop Solution) to each well.
13. Determine the optical density within 30 minutes, using microplate **reader** set to 450 nm corrected with 570 nm or 630 nm.

**TYPICAL DATA**

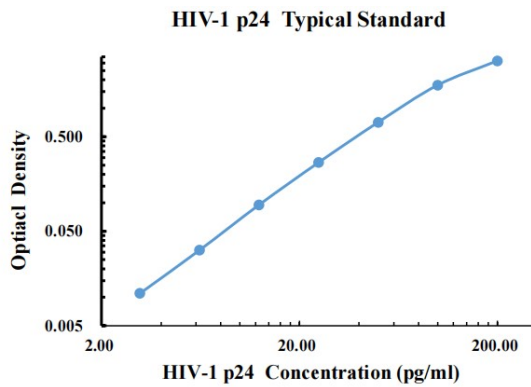

pg/ml	O.D.		Average	Corrected
0.00	0.0190	0.0200	0.0195	
3.13	0.0300	0.0310	0.0305	0.0110
6.25	0.0520	0.0500	0.0510	0.0315
12.50	0.1150	0.1140	0.1145	0.0950
25.00	0.2920	0.2800	0.2860	0.2665
50.00	0.7590	0.7000	0.7295	0.7100
100.00	1.8220	1.7240	1.7730	1.7535
200.00	3.2790	3.0910	3.1850	3.1655

**RECOVERY**

The spike recovery was evaluated by spiking 3 levels of human HIV-1 p24 into health human serum sample. The un-spiked serum was used as blank in this experiment.

The recovery ranged from 94% to 105% with an overall mean recovery of 98%.

**LINEARITY**

To assess the linearity of the assay, five samples were spiked with high concentration of HIV-1 p24 in human serum and diluted with Sample Diluent to produce samples with values within the dynamic range of the assay.

The linearity ranged from 95% to 104% with an overall mean recovery of 99%.

**SAMPLE DILUENT**

If the concentration of the sample is too high, the sample can be diluted with **Sample Diluent**.

**DISCLAIMER AND VERSION**

For research use only. Not for use in clinical diagnostic procedures.

Copyright 2024 AREX Biosciences Ltd. Version.240820

**SENSITIVITY**

The minimum detectable dose (MDD) of human HIV-1 p24 is typically less than 0.64 pg/ml (50 µl of sample volume) or 2.62 pg/ml (10 µl of sample volume).

The MDD was determined by adding two standard deviations to the mean optical density value of ten zero standard replicates and calculating the corresponding concentration.

**PRECISION**

- Intra-assay Precision (Precision within an assay) Three samples of known concentration were tested twenty times on one plate to assess intra-assay precision.
- Inter-assay Precision (Precision between assays).

	Intra-assay Precision			Inter-assay Precision		
	S1	S2	S3	S1	S2	S3
Sample Number	22	22	22	6	6	6
Average (pg/ml)	4.6	18.4	58.9	4.2	19.5	59.7
Standard deviation	0.12	0.59	4.1	0.21	0.72	3.8
Coefficient of variation (%)	2.6	3.2	6.9	5.0	3.7	6.4